


Our flagship in
urinalysis




- Fully automated urine analyzer
- User-friendly touch screen
- Throughput of 250 tests per hour
- Interfaces for connection to laboratory network
- Integrated barcode reader

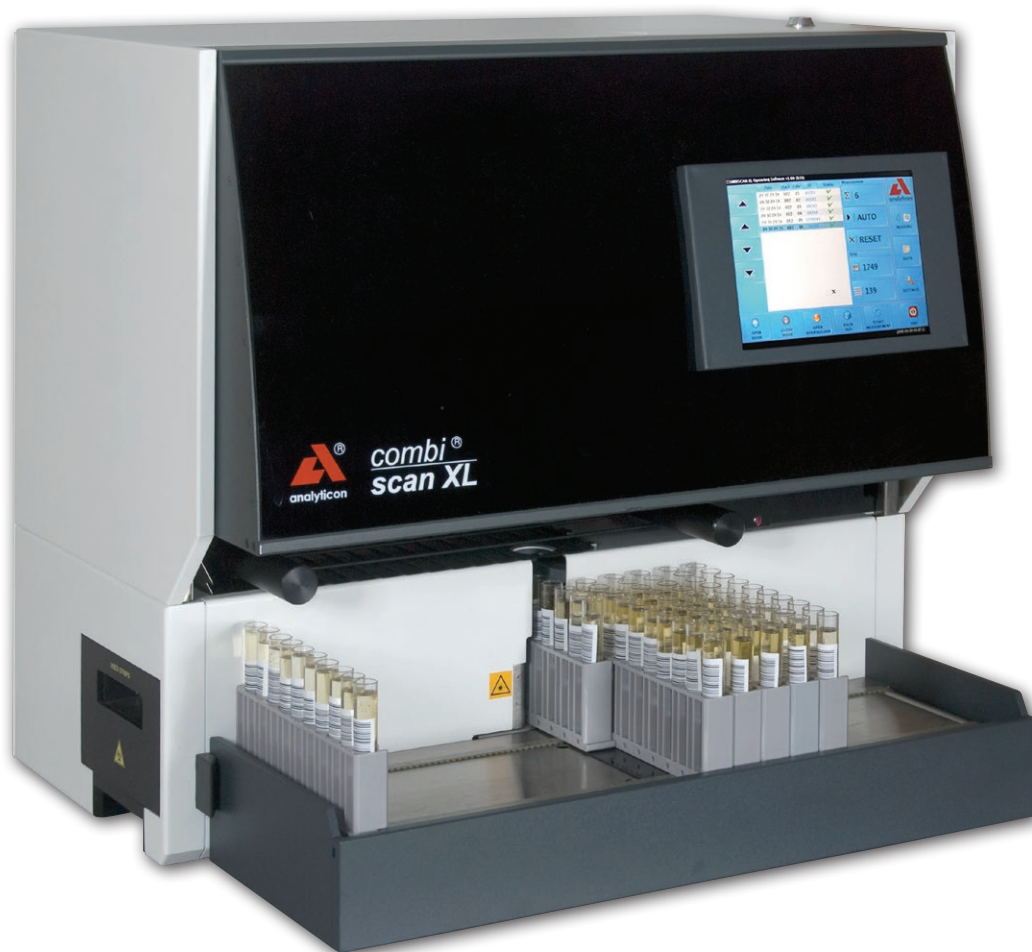

Internal test strip container for 150 test strips



Waste bin for hygienic disposal of used test strips

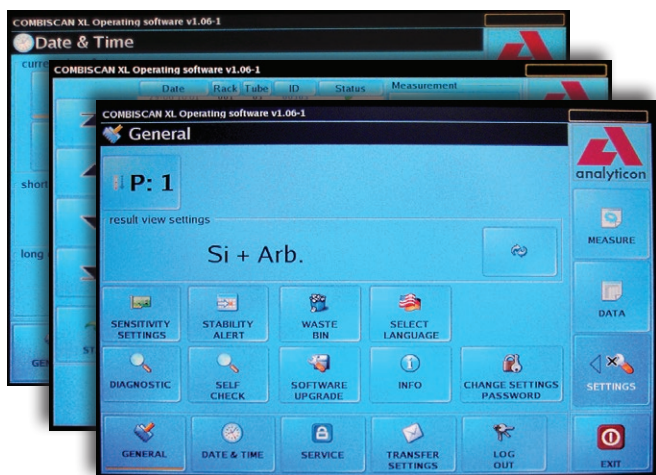

User-friendly high resolution touch screen

Features & Functions:

- Fully automated test strip and sample handling
- Big touch screen for convenient operation
- Integrated barcode reader
- Data entry via keyboard
- Loading of up to 100 samples possible
- Urine colour compensation
- Interfacing with laboratory network system possible


Technical data:

- Type: CCD Colour Image Sensor
- Throughput: 250 strips/hour
- Display: 170 x 130 mm touch screen
- Memory: Last 9999 results
- Interfaces: USB, RS 232, serial interface, monitor, LAN

Optionals:

- Keyboard, monitor, mouse
- Extension by sediment analysis system possible (under development)

(All modifications in view of technical progress are reserved)

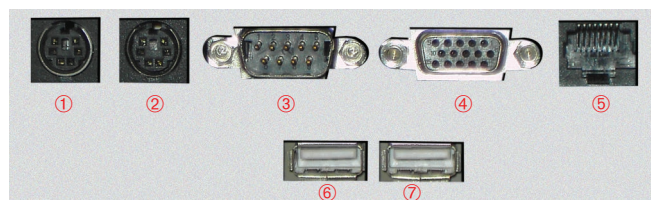

Reverse side connectors: ① = keyboard, ② = mouse, ③ = serial interface, ④ = monitor, ⑤ = LAN, ⑥ + ⑦ = USB



Urine Strip Analyzer

Combi Scan® XL

Art.-#

A93 007

Urine Test Strips

Combi-Screen® 11SYS

Combi-Screen® 11SYS PLUS

Pack size

100 / 150 tests

100 / 150 tests

Art.-#

93 100 / 93 150

94 100 / 94 150

Parameters:

Bilirubin, Urobilinogen, Ketones, Ascorbic Acid, Glucose, Protein, Blood, pH, Nitrite, Leukocytes, Specific Gravity

Combi-Screen® and Combi Scan® are registered trademarks of Analyticon Biotechnologies AG.

Distributor information

**Analyticon
Biotechnologies AG**

Am Muehlenberg 10
35104 Lichtenfels - Germany

Fon +49 (0) 6454 7991-0
Fax +49 (0) 6454 7991-71

info@analyticon-diagnostics.com
www.analyticon-diagnostics.com